

The Use of CD45/SSC Dot Plots in the Classification of Acute Leukemias

Mohamed Brahimi^{1,2*}, Djamel Saidi³, Hadj Touhami³ and Mohamed Amine Bekadja²

¹Department of Hemobiology, EHU "1er Novembre 1954" Oran, Algeria

²Department of Hematology and Cell Therapy, EHU "1er Novembre 1954" Oran, Algeria

³Department of Hematology, CHU Oran, Algeria

The anti-CD45 antibody is the most used reagent for immunophotyping acute leukemias. The analysis of CD45/SSC dot plot is mandatory when analyzing leukemic cells because this graph allows cytometristes to isolate blast cells from other hematopoietic cell types.

CD45 also called "Leukocytes Common Antigen" is expressed on mature and immature hematopoietic cells with the exclusion of platelets, erythrocytes, elements belonging to the erythroid compartment, and highly immature hematopoietic precursors. It is well known that lymphocytes are brightly stained with CD45 and gives a weak SSC

signal, monocytes got less CD45 on their membrane than lymphocytes with a slightly higher SSC signal while granulocytes are weakly stained with CD45 and gives a very high SSC compared to the other nucleated cells.

On a CD45/SSC bi-dimensional plot blast cells are located in the so-called "Bermude Area"; in the middle of the other described cell types. According to Haycocks, CD45/SSC cytogram allowed the flow cytometristes to classify immature leukemic cells into lymphoblasts and myeloblasts in 94, 1% of cases. In their paper the authors described the blasts of AML as a localized round/oval cluster with slightly higher SSC signal than lymphocytes (Figure 1) while the blasts of ALL as clusters distributed longitudinally along the CD45 axis with SSC signal similar to lymphocytes (Figure 2). Monoblastic population gives a very heterogeneous pattern which overlaps with the monocytic area [1].

Since different nucleated cells in a lysed whole blood and bone marrow can be quantified in a CD45/SSC plot, it is possible to make a correlation between the cytometric pattern and the FAB subtypes in AML. Figure 3 shows typical CD45/SSC patterns for in each AML FAB subtype according to Saidi et al. [2], but it should not be forgotten that on the promyelocytes of AML-M3v, CD45 is expressed at a higher intensity than on the promyelocytes of the classic form. In contrast, it is also possible that the blasts of sporadic cases of AML do not express the antigen; in such cases the immunological gate must be modified depending on the blast cell phenotype [3].

On T-lymphoblasts, CD45 is usually expressed at a lower intensity than on normal T lymphocytes, but at a somewhat higher intensity than on the cells of neoplastic diseases of B cell precursors. Consequently, the cluster of T leukemic lymphoblasts can partially overlap with the cluster of residual normal lymphocytes [3].

On the blasts of B-ALL, the intensity of CD45 expression is lower than on mature lymphocytes, and similar to the normal B precursors. But the intensity of CD45 on lymphoblasts is more heterogynous than on myeloblast that is why they are spread along the CD45 axis. However some special features of CD45/SSC in B-ALL must be kept in mind such as: 1) In cases with hyperdiploidy or with translocation t (12;21), CD45 can be missing or expressed at an intensity similar to erythroblasts. 2) In some cases with translocation t (4;11), the intensity of CD45 expression can be similar to mature lymphocytes. 3) Finally, in some cases with translocation t(9;22), a high expression of CD45 together with a high expression of CD19, CD22, CD34 and HLA-DR correlates with the presence of trisomy 8 [3].

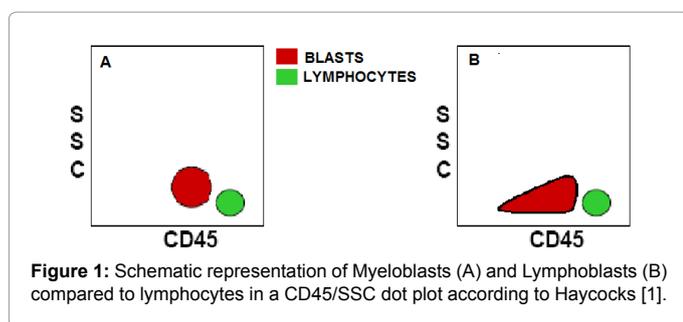


Figure 1: Schematic representation of Myeloblasts (A) and Lymphoblasts (B) compared to lymphocytes in a CD45/SSC dot plot according to Haycocks [1].

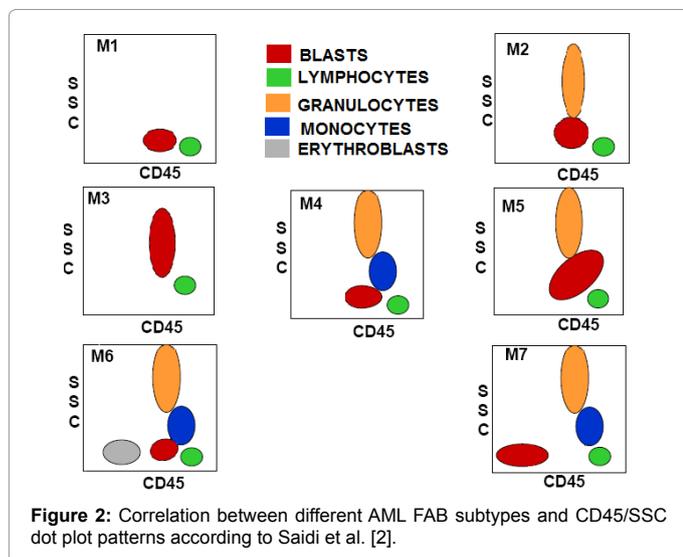


Figure 2: Correlation between different AML FAB subtypes and CD45/SSC dot plot patterns according to Saidi et al. [2].

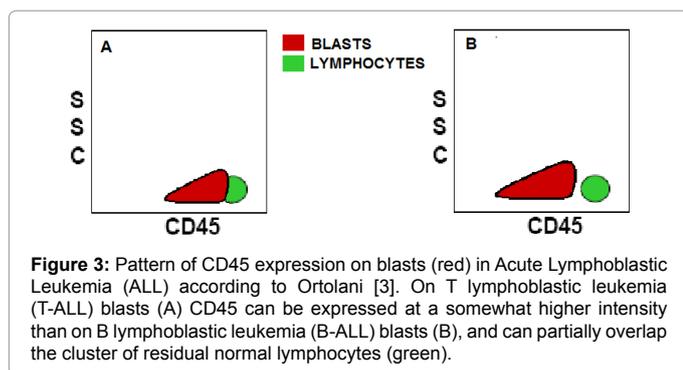


Figure 3: Pattern of CD45 expression on blasts (red) in Acute Lymphoblastic Leukemia (ALL) according to Ortolani [3]. On T lymphoblastic leukemia (T-ALL) blasts (A) CD45 can be expressed at a somewhat higher intensity than on B lymphoblastic leukemia (B-ALL) blasts (B), and can partially overlap the cluster of residual normal lymphocytes (green).

***Corresponding author:** Mohamed Brahimi, Department of Hematology and Cell Therapy, Etablissement Hospitalier et Universitaire (EHU) "1er Novembre 1954" Oran, Algeria, Tel: 00213 771 485739; Fax: 00213 41 42 16 36; E-mail: bmw73dz@yahoo.fr

Received January 19, 2014; Accepted January 23, 2014; Published January 25, 2014

Citation: Brahimi M, Saidi D, Touhami H, Bekadja MA (2014) The Use of CD45/SSC Dot Plots in the Classification of Acute Leukemias. J Hematol Thromb Dis 2: e107 doi: [10.4172/2329-8790.1000e107](https://doi.org/10.4172/2329-8790.1000e107)

Copyright: © 2014 Brahimi M, et al. This is an open-access article distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited.

References

1. Haycocks NG, Lawrence L, Cain JW, Zhao XF (2011) Optimizing antibody panels for efficient and cost-effective flow cytometric diagnosis of acute leukemia. *Cytometry B Clin Cytom* 80: 221-229.
2. Saidi D, Brahimi M, Mekkous-Touhami F, Zatta L, Mrabet R, et al. (2007) L'immunophénotypage par cytométrie en flux des leucémies aiguës non lymphoblastiques. 4^{ème} Congrès National d'Hématologie Alger.
3. Ortolani C (2007) *Flow Cytometry of Hematological Malignancies*. (1stedn), John Wiley & Sons Ltd., UK.

Citation: Brahimi M, Saidi D, Touhami H, Bekadja MA (2014) The Use of CD45/SSC Dot Plots in the Classification of Acute Leukemias. J Hematol Thromb Dis 2: e107 doi: [10.4172/2329-8790.1000e107](https://doi.org/10.4172/2329-8790.1000e107)

Submit your next manuscript and get advantages of OMICS Group submissions

Unique features:

- User friendly/feasible website-translation of your paper to 50 world's leading languages
- Audio Version of published paper
- Digital articles to share and explore

Special features:

- 300 Open Access Journals
- 25,000 editorial team
- 21 days rapid review process
- Quality and quick editorial, review and publication processing
- Indexing at PubMed (partial), Scopus, EBSCO, Index Copernicus and Google Scholar etc
- Sharing Option: Social Networking Enabled
- Authors, Reviewers and Editors rewarded with online Scientific Credits
- Better discount for your subsequent articles

Submit your manuscript at: <http://www.omicsonline.org/submission/>

